



fmGenie Neck Loop

1a


fmGenie

For a hearing aid to be able to 'listen' to the fmGenie, you need to use either a neck loop or direct input connection.

An inductive neck loop can be used with any hearing aid which has a 'T' or Telecoil facility.

A neck loop can be worn under clothing and whilst the sound quality may not be as good as using a direct input, the use of an inductive neck loop is sometimes preferred on the cosmetic grounds of not wanting to have wires visible – and some hearing aids do not have a direct input facility anyhow.

With the majority of hearing aids when you switch to the 'T' position you will mute the hearing aid microphone. To overcome this loss of own voice you will also need a microphone pack for the fmGenie receiver – see page 15.

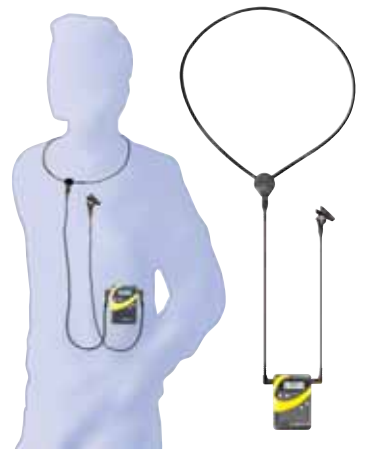

Multiturn neck loop


Multiturn plastic neck loop

This multiturn plastic loop has a soft flexible feel and lies comfortably around the neck whether under or over your clothes.

This neck loop has a 2 pin N.P. socket and requires the use of the patch cord shown below.

Plastic neck loops are approx. 225mm diameter. Thickness of loop: 5.5mm. Weight: 35g.

FMG441 Multiturn plastic neck loop £ 24.36 each

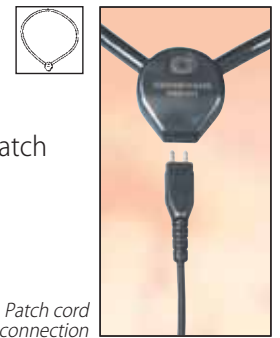

Patch cord connection

Patch cords for fmGenie multiturn plastic neck loop

Black with a straight 2 pin N.P. plug for plugging into the neck loop and a right angled 3.5mm plug to connect to the receiver. These neck loop patch cords are available in a range of lengths and can also be used with bone conductors – see page 25.

- FMG42200 200mm standard patch cord £ 8.76 each
- FMG42400 400mm standard patch cord £ 8.96 each
- FMG42600 600mm standard patch cord £ 9.17 each
- FMG42800 800mm standard patch cord £ 9.37 each



What the ends look like!



Please note that this neck loop is not recommended for use in normal stereo headphone sockets as it will only pick up one channel. A personal inductive stereo neck loop, which looks similar but which offers excellent sound quality, can be found in Section 11.



Listening to an fmGenie with headphones

Using headphones with an fmGenie receiver allows people who have listening and concentration difficulties to use an fmGenie system to benefit from sound reinforcement.

fmGenie headphones

Personal stereo style headphones for users requiring help with concentration in lectures or for those with a mild hearing loss.

FMG7420 £ 5.56 each



fmGenie test headphones

A lightweight pair of headphones for routine listening tests, not particularly robust but comfortable to wear with a pleasant sound; very good value for money!

FMG7410 £ 2.14 each

